




Maska
Face masks for anaesthesia

Maska face masks are ideal for use with all types of manual and automatic resuscitators or for anaesthesia. The structure is anatomically and ergonomically modelled to guarantee a perfect fit on the patient's face and an optimum grip for the operator.

The masks are antistatic and have an inflatable cushion to ensure a perfect fit on the patient's face. The masks are made of rubber and they are supplied with a stainless steel ring for nape fixation.



- Maska
 - Maska
 - Maska
 - Maska
 - Maska
 - Maska
- size 0 Code RM20600B
 - size 1 Code RM20601B
 - size 2 Code RM20602B
 - size 3 Code RM20603B
 - size 4 Code RM20604B
 - size 5 Code RM20605B

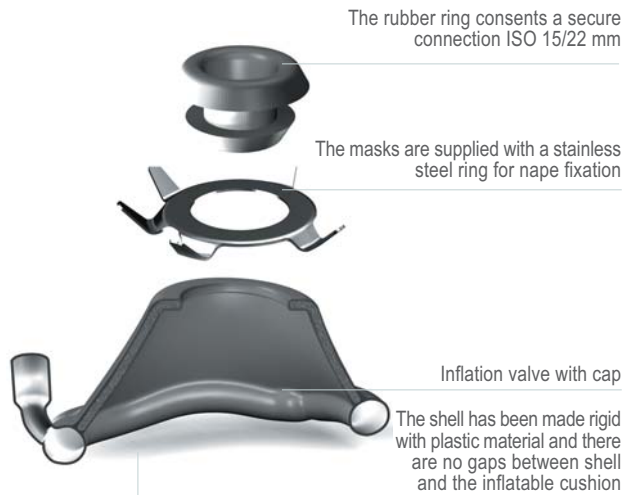



Spencer Mask
Face masks with inflatable cushion

Spencer transparent face masks are ideal for use with all types of manual or automatic resuscitators. The structure is modelled anatomically and ergonomically to guarantee a perfect fit to the patient's face and an optimum grip for the operator.

Inflatable cushion

The inflatable cushion is made of silicone and does not deflate: the mask is always ready to be used. The pressure inside the cushion can be modified, in order to adapt and guarantee a good fit to the patient's face. The use of silicone has enabled us to present the first latex free inflatable mask.



- Spencer Mask black poly-sulphonate
 - size 00 Code RM20806A
 - size 0 Code RM20800A
 - size 0 low Code RM20801A
 - size 2 Code RM20802A
 - size 4 Code RM20804A
 - size 5 Code RM20805A
- Spencer Mask black poly-carbonate
 - size 00 Code RM20706A
 - size 0 Code RM20700A
 - size 0 low Code RM20701A
 - size 2 Code RM20702A
 - size 4 Code RM20704A
 - size 5 Code RM20705A

The silicone ring permits a secure ISO connection 15/22 mm



The masks can also be supplied with a stainless steel ring for nape fixation

The shell is made of poly-sulphonate or poly-carbonate and the transparent shell permits the inspection of the patient's spontaneous breathing



Assemblage is simple thanks to buttons/knobs on the shell



The valve is connected to the cushion with a specially designed ring



The inflatable silicone cushion is modelled anatomically and it has a soft surface to increase the patient's comfort



Special part of the inflatable cushion



**B-life
Silicone face masks**

B-life face masks are ideal for use in anaesthesia and in assisted and controlled ventilation. They are autoclavable and characterised by maximum adhesion. The connecting ring is soft to guarantee easy attachment and to avoid disconnections.


The shape of the transparent shell makes them comfortable for the operator and also on the patient's face. It is supplied with a standard ISO connection 15/22 mm and minimal dead space.




- B-life Silicone mask size 0 Code RM20400B
- B-life Silicone mask size 1 Code RM20401B
- B-life Silicone mask size 2 Code RM20402B
- B-life Silicone mask size 3 Code RM20403B
- B-life Silicone mask size 4 Code RM20404B
- B-life Silicone mask size 5 Code RM20405B

Available in rubber


- Eco B-life Rubber mask size 0 Code RM20504B
- Eco B-life Rubber mask size 1 Code RM20503B
- Eco B-life Rubber mask size 2 Code RM20502B
- Eco B-life Rubber mask size 3 Code RM20510B
- Eco B-life Rubber mask size 4 Code RM20501B
- Eco B-life Rubber mask size 5 Code RM20505B



The silicone ring allows secure ISO connection 15/22 mm
In rubber for Eco B-life model



The shell is made of transparent poly-carbonate, in order to monitor the patient's spontaneous respiration



The silicone cushion of the adolescent and adult version is modelled to improve the fit to the patient's face and thanks to incorporated buttons, assemblage is very quick after sterilisation
In rubber for Eco B-life model



Air Cuffed Mask
PVC disposable face masks

The air cuffed masks can be used with manual and automatic resuscitators and ventilators. The design and the soft pneumatic ring guarantee maximum adaptability to the shape of the face. The use of economical materials (the shell and the cushion are made of PVC) promotes the idea of using the mask for a single patient, avoiding cross infection and high sterilisation costs.

Complete with plastic ring

They are available in 8 sizes, with or without inflation valve. Standard ISO connection 15/22 mm.

Dead space:

- adult sizes: 150 ml
- paediatric sizes: 80 ml
- neonatal sizes: 40 ml



- | | |
|---|----------------------|
| - Air Cuffed Mask | size 0 Code RM10880A |
| - Air Cuffed Mask | size 1 Code RM10881A |
| - Air Cuffed Mask | size 2 Code RM10882A |
| - Air Cuffed Mask | size 4 Code RM10884A |
| - Air Cuffed Mask | size 5 Code RM10885A |
| - Air Cuffed Mask colour coded with valve | size 0 Code RM10900A |
| - Air Cuffed Mask colour coded with valve | size 1 Code RM10901A |
| - Air Cuffed Mask colour coded with valve | size 2 Code RM10902A |
| - Air Cuffed Mask colour coded with valve | size 3 Code RM10903A |
| - Air Cuffed Mask colour coded with valve | size 4 Code RM10904A |
| - Air Cuffed Mask colour coded with valve | size 5 Code RM10905A |



Plastic cross for the nape fixation, colour-coded for every size

The PVC ring allows secure ISO connection 15/22 mm



The shell is made of transparent PVC, in order to monitor the patient's spontaneous respiration

The soft inflated cushion, which is joined to the shell assures an optimum fit to the patient's face



Pro Valve
One way valve for CPR

It allows you to practice "mouth to mouth" resuscitation. Thanks to the non-rebreathing valve welded with ultrasonic weldings on a transparent PVC sheet, it avoids infection transmission. The water-repellant filter avoids contact with patients' saliva or blood.
Made of non-toxic transparent material, easy to use and suitable for adults and children. It can also be used on manikins .

Pro Valve

Code TA07001C



Pro Valve Key 1
Nylon key-ring

Soft and handy key-ring, made of tear-proof nylon Ripstop with Star of Life symbol, it always allows you to bring along the Pro-Valve.

Pro Valve Key 1
Pro Valve Key 1 with Pro Valve

Code TA07012C
Code TA07011C



Key-life bag 2
Nylon Key-ring

Key-ring in tearproof material. An extremely small space always allows you to bring along this vital device.

Key-life bag 2 blue
Key-life bag 2 red
Key-life bag 2 orange
Key-life bag 2 yellow

Code CT06070B
Code CT06071B
Code CT06072B
Code CT06073B



Flex Mask
One way mouth to mouth mask

Spencer Flex Mask comprises a one-way valve, which is connected to a resuscitation mask, made in a latex free material. It prevents mouth to mouth contact. Flex Mask can be used any time a situation arises, in which a patient in respiratory arrest has to be resuscitated. The air flow can be enhanced with oxygen by removing the cap on the mask and by inserting the oxygen supply tube. The transparent mask allows the operator to immediately see any strange objects that may be regurgitated by the patient during resuscitation. The softness of the cushion guarantees the perfect adhesion to the patient's face. Flex Mask is a must have for every first aid kit and above all... pocket! Flex Mask, as well as being a life saving instrument, can also be used on training mannequins.

Flex Mask with orange PVC box with logo and printed instructions
Flex Mask without valve
Spare valve
Disposable filter

Code TA09019A
Code TA09010A
Code TA09030A
Code TA09040A



EVA Mask Case
Soft and elegant container in EVA equipped with loop with automatic closure for belts

Flex Mask with blue EVA box

Code TA09050A

PVC Mask Case
The mask can be stored in an unbreakable, polypropylene, tapered, light and compact case, complete with snap closure.

Flex Mask with orange PVC box with logo and printed instructions
Flex Mask with red PVC box with logo and printed instructions
Flex Mask with yellow PVC box with logo and printed instructions
Flex Mask with green PVC box with logo and printed instructions
Flex Mask with blue PVC box with logo and printed instructions

Code TA09019A
Code TA09033A
Code TA09034A
Code TA09035A
Code TA09036A

Nylon Mask Case
Flex Mask can be stored in a case made of unbreakable nylon with a belt loop. The version with latex gloves and antimicrobial serviette contributes to optimising the operators protection.

It is supplied with:
- two serviettes for medication
- one pair of gloves

Flex Mask with yellow/black material mask case

Code TA09037A